




Distillation of Cutback Asphaltic Products

Shield
Internal diameter 117 mm
Total height 117 mm
Mica windows 45 x 45 mm
Rim external diameter 148 mm

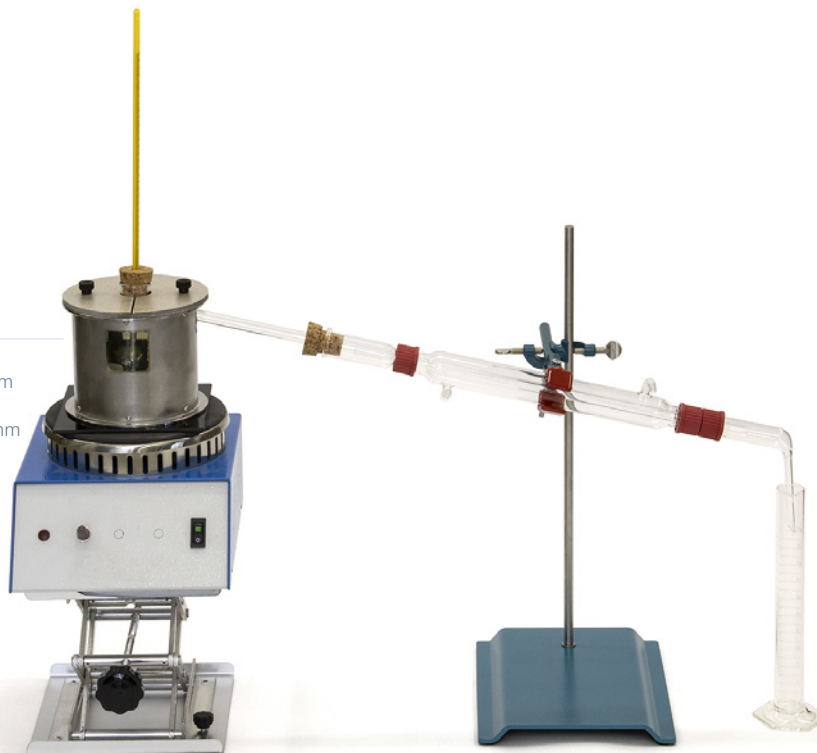

ASTM D402

Distillation of Cutback Asphaltic (Bituminous) Products

LT/CB-106000/M

Manual instrument composed by:

- Electric heater 1000 Watt with main switch and power regulator, supported on a manual height adjustable platform
- Shield made in steel lined with 3 mm fire proof insulation and fitted with transparent mica windows, cover in two halves made of fire resistant material
- Distillation flask, 500 ml side-arm with stopper for thermometer
- Nozzle extensor made in glass
- Water condenser made in glass with 8 mm joints
- Metal base with rod, clamp and clamp for glassware
- Adapter made in glass with reinforced top for avoiding splash
- Receiver cylinder graduated to 100 ml, div. 1 ml

Power supply

- 220 or 115 Vac 50/60 Hz

Accessories

- T-AS8C: Thermometer ASTM 8C - IP 6C

Spare Parts

- LAB-101-050: distillation flask, 500 ml side-ACM with stopper for thermometer
- LAB-101-051: shield made in steel lined with 3 mm fire proof insulation and fitted with transparent mica windows, cover in two halves made of fire resistant material
- LAB-101-052: cover in two halves
- LAB-101-053: set of stoppers
- LAB-101-054: receiver cylinder 100 ml, div. 1 ml
- LAB-101-055: water jacketed condenser made in glass
- LAB-101-056: adapter made in glass
- LAB-101-057: nozzle extensor 450 mm